Organic Hemp CO2 Extract

Sunshine Trading | PHYTO Plus Zwanebloemlaan 222 1087GD Amsterdam NH The Netherlands

Phone: +31(0)20 770 97 91 support@phytopluscbd.com www.phytopluscbd.com Responsible Supervisor: Martin V.

Sample Batch #3404

Date samples received: 02 December 2021
Date analysis began: 02 December 2021
Date analysis began: 04 December 2021

Date sample report produced: 06 December 2021

ID Number when available:

Sample Mass 1 g






SKAL:100364, ISO 22000 certified, Organic certified: NL-BIO-01, HACCP certified; GMP certified

Phyto Plus 6.14% Total CBG: Cannabinoid Profile


Component	Mass (%)	Amount (mg/g)
CBD	0.91	9.10
CBDA	0.07	0.70
CBDV	<loq< td=""><td><loq< td=""></loq<></td></loq<>	<loq< td=""></loq<>
CBG	6.14	61.40
CBGA	0.05	0.50
CBN	0.07	0.70
CBC	0.02	0.20
Δ9ΤΗС	<loq< td=""><td><loq< td=""></loq<></td></loq<>	<loq< td=""></loq<>
THCA	<loq< td=""><td><loq< td=""></loq<></td></loq<>	<loq< td=""></loq<>
Δ8THC	<loq< td=""><td><loq< td=""></loq<></td></loq<>	<loq< td=""></loq<>
Total CBG	6.14	61.40

LOQ - Limit of Quantitation (LOQ=0.03%, LOD=0.006%)

Method: HPLC-UV

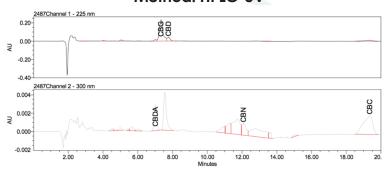

Cannabinoids as Percent of Total Mass

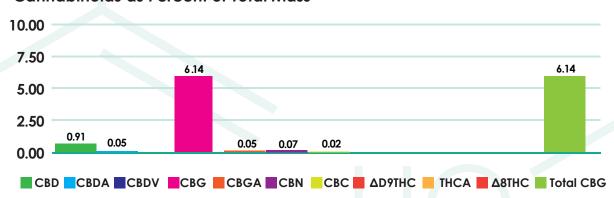

Organic Hemp CO2 Extract

Terpenoid Profile

Sunshine Trading | PHYTO Plus Zwanebloemlaan 222 1087GD Amsterdam NH The Netherlands

Phone: +31(0)20 770 97 91 support@phytopluscbd.com www.phytopluscbd.com Responsible Supervisor: Martin V.
Sample Batch #3404
Date samples received: 02 December 2021
Date analysis began: 02 December 2021
Date sample report produced: 06 December 2021
ID Number when available:




SKAL:100364, ISO 22000 certified, Organic certified: NL-BIO-01, HACCP certified; GMP certified

Sample Mass 1 g


Phyto Plus 6.14% Total CBG: Terpenoid Profile

Component	Amount (%)
β-Caryophyllene	0.07
a-Humulene	<0.01
Caryophyllene oxide	<0.01
Myrcene	0.31
a-Pinene	0.69
Terpinolene	<0.01
Humulene epoxide II	<0.01
Limonene	0.83
β-Pinene	<0.01
E-β-Ocimene	<0.01
Sabinene	<0.01
Linalool	0.13

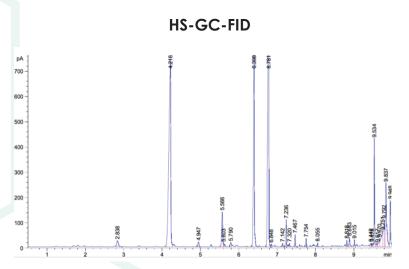

Terpenoid Distribution



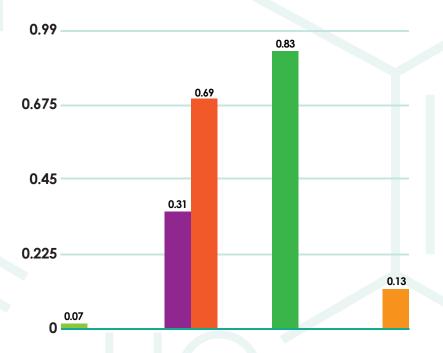

Organic Hemp CO2 Extract

Sunshine Trading | PHYTO Plus Zwanebloemlaan 222 1087GD Amsterdam NH The Netherlands

Phone: +31(0)20 770 97 91 support@phytopluscbd.com www.phytopluscbd.com Responsible Supervisor: Martin V.
Sample Batch #3404

Date samples received: 02 December 2021
Date analysis began: 02 December 2021
Date sample report produced: 06 December 2021

ID Number when available:

Sample Mass 1 g




SKAL:100364, ISO 22000 certified, Organic certified: NL-BIO-01, HACCP certified; GMP certified

Phyto Plus 6.14% Total CBG: Microbial Profile

Component	Amount (mg/g)	Results
Listeria m.	1 g	ND
Escherichia c.	1 g	ND
Salmonella	25 g	ND
Yeast	1 g	ND
Mould	1 g	ND

*ND - Not detected

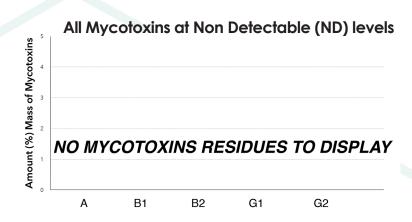

Nutrition Facts

Component	%
Moisture	2.15
Protein	0.31
Total fat	97.25
Total Carbohydrates	ND
Dietary Fibers	ND
Sugars	ND
Ash	ND

^{*}ND - Not detected

Organic Hemp CO2 Extract

Sunshine Trading | PHYTO Plus Zwanebloemlaan 222 Sample Batch #3404 1087GD Amsterdam NH The Netherlands


Phone: +31(0)20 770 97 91 support@phytopluscbd.com www.phytopluscbd.com

Responsible Supervisor: Martin V.

Date samples received: 02 December 2021 Date analysis began: 02 December 2021 Date sample report produced: 06 December 2021

ID Number when available:

Sample Mass 1 g

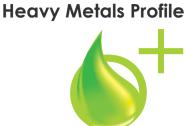

SKAL:100364, ISO 22000 certified, Organic certified: NL-BIO-01, HACCP certified; GMP certified

Phyto Plus 6.14% Total CBG: Heavy Metals Profile

Component	Mass (%)	Amount (ppm)	Limit (ppm)
Arsenic (As ₂ O ₃)	ND	< 0.1	< 0.1
Cadmium (Cd)	ND	< 0.1	< 0.1
Lead (Pb)	ND	< 0.1	< 0.1
Mercury (Hg)	ND	< 0.1	< 0.1
Chromium (Cr)	ND	< 1	< 1
Tin (Sn)	ND	< 10	< 10

ND - Not detected, **Codex STAN 193-1995, GB 2762, EC No. 1881/2006, FDA

All Heavy Metals at Non Detectable (ND) levels



Conclusions:

No heavy metal residues detected. No flammable residues detected. No chemical residues detected.

Organic Hemp CO2 Extract

Pesticide Analysis

Sunshine Trading | PHYTO Plus Zwanebloemlaan 222 1087GD Amsterdam NH The Netherlands

Phone: +31(0)20 770 97 91 support@phytopluscbd.com www.phytopluscbd.com Responsible Supervisor: Martin V.
Sample Batch #3404
Date samples received: 02 December 2021
Date analysis began: 02 December 2021
Date sample report produced: 06 December 2021

ID Number when available:

Sample Mass 1 g



SKAL:100364, ISO 22000 certified, Organic certified: NL-BIO-01, HACCP certified; GMP certified

Pesticide Analysis: Our tests looked for residue of nearly 300 known pesticides finding no evidence of any over detectable limits.

The Lab tests our products thoroughly. Nearly 300 of the below pesticides concentrations were measured and we are proud to say that all tests measured below our detectable limits. Most tests have a threshold of 0.01 mg/k, while only a handful of tests have a threshold value of <0.05 mg/kg. Not a single test of PHYTO Plus products went over detectable threshold limits.

PESTICIDES MEASURED

Acrinathrin Azoxystrobin Biphenhin Bitertanol Biphenyl Bromopropylate Bromuconazole Bupirimate Cadusafos Captafol Captan Chlorphenson Chlorfenapyr Chlorfenvinphos Chlorothalonil Chlorprophame 3,5-Dichloraniline Chlorpyriphos Chlorpyrifos-methyl Chlorthal-dimethyl Cyfluthrin Cypermethrin Cyproconazole Cyprodinil Clomazone o,p-DDE P,P-DDE o,p-DDD P,P-DDD o,p-DDT p,p-DDT Deltamethri Diazinon Diclofop-methyl Dieldrin Dichlobenil Dichlofluanid Dichloryos Dicloran Dicofol Dicrotophos Diethofencarb Diflubenzuron Dimetachlor Diniconazole Dodemorph Diphenylamine Alpha-Endosulfan Beta-Endosulfan Endosulfan-sulphate Ethion Etofumesate Ethoprophos Ehtoxyquin Etoxazole Etridiazole Etrimphos Famoxadone Fenarimol Fenazaquin Fenchlorphos Fenhexamid Fenihothion Fenpropidin Fenpropimorph Fenvalerate Formothion Fipronil Fipronil-sulfone Fludioxonil Flusilazole Flutriafol Folpet Fuberidazole Furathiocarb Hexaconazole HCB Alpha-HCH Beta-HCH Delta-HCH Heptachlor Heptachlor-epoxidceis Heptachlorepoxidtreans Iprodione Iprovalicarb Lambda- cyhalothrin Lindane Mecarbam Metalaxv Metazachlor Methidathion Metribuzin Mevinphos Myclobutanil Nuarimol Orthophenylphenol Oxadixyl Paclobutrazol Parathion Parathion-methyl Paraoxon-methyl Paraoxon-ethyl Penconazole Pendimethaline Permethrin Phenthoate Phorate Procymidone Profenofos Propiconazole Propyzamide Pyrazophos Pyrethrins Pyridaben Pyrimethanil Pyriproxyfen Quinoxifen Quitozene Pentachloraniline Phosphamidon Pyrifenox Prometryn Propanil Propoxur Proquinazid Prothiofos Simazine Spiroxamine T au-fluvalinate T ebuconazole T ebufenpyrad T ecnazene T efluthrin T erbuthylazine T etraconazole T etradifon T etramethrine T olclofos-methyl T olylfluanid Transfluthrin Triadimephon Triadimenol Trialate Trifloxystrobin Triflumizole Vinclozolin DDT isomersum Heptachlor (heptachloarnd heptachloer poxidsum) Trifluraline Chlorobenzilate 3-Chloraniline Abamectin (AvermectinBla and AvermectinBlb sum) Acetamiprid Aldicarb Aldikarbsulphone Aldicarbsulphoxide Azinphos-ethyl Azinphos-methyl Benalaxyl Benfuracarb Boscalid Buprofezin Carbaryl Carbendazim Carbofuran 3-hydroksicarbofuran Carbosulfan Chloridazon Cymoxanil Clofentezin Clothianidin Demeton-S-methyl Demeton-S-methyslulfoxid Diafenthiuron Difenoconazole Dimethoate Dimethomorph Diuron EPN Epoxiconazole Ethirimol Etofenprox Fenamidone Fenbuconazole Fenbutatinoxid Fenoxycarb Fenpyroximate Fenpropathrin Fensulfothion Fenthion Fenthionsulphone Fenthionsulphoxide Fluazinam Flufenoxuron Fluquinconazole Fonofos Formetanate Fosthiazate Hexythiazox Imazalil Imidacloprid Indoxacarb Isofenphos Methacrifos Isofenphosmethyl Krezoxim-methyl Linuron Lufenuron Malaoxon Malathion Mepanipirim Mepronil Metamitron Metconazole Methamidophos Methiocarb Methiocarbsulphone Methiocarbsulfoxide Methomyl Methoxyfenozide Metobromuron Monocrotophos Monolinuron Omethoate Oxamyl Pencycuron Phenmedipham Phosalone Phosmet Phosmeot xon Phoxim Pymetrozine Piperonylbutoxide Pyraclostrobin Pyridaphenthion Pyridate Pyrifenox Pirimicarb Pirimicarbdesmethyl Pirimiphos-methyl Primisulfuron-methyl Prochloraz Propamocarb Propargite Prothioconazole Prothioconazole-desthio Quinalphos SpinosynA SpinosynD Sulfotep T ebufenozide T eflubenzuron Thiabendazole Thiacloprid Thiamethoxam Thiodicar Thiophanatemethyl Tralkoxydim Triazophos Trichlorfon Triflumuron Triforine Triticonazole Zoxamide Acephate Amitraz Fenamiphos Fenamiphosulphone Fenamiphosulfoxid Nitempiram Fenthionoxonsulphone Fenthionoxonsulfoxid Kumapho Piriphenox Mehibuzine DEET

Our laboratory analysis is standardized after following protocols:

ISO 6579:2003 ISO 11290-1:2003 ISO 16649-2:2002 ISO 21527-2:2008

Note on Cannabinoid Testing:

All cannabinoids in their acid forms (ending in "-A") are convertible to their non-acid forms via a decarboxylation process (heating). The components lose mass through this process. To find the total theoretical active cannabinoids, one multiplies the acid forms by 87.7%. For example, THC-A can be converted to active THC using the formula: THC-A x 0.877 = THC. In this case, the Max THC for the sample is: Max THC = (THC-A x 0.877) + THC. This method has been validated according to the principles of the International Conference on Harmonisation.

Chromatographic Analysis:

Analysis of cannabinoids content was performed using Waters 2695 (Milford, MA, USA) separation module equipped with auto injector, sample cooler, vacuum degasser and column heater units. Separation of all cannabinoids was accomplished on YMC PRO C18 (150 x 4 mm I.D., S-3µm) RP column coupled with C18 precolumn maintained at 30 °C by a CTO-20AC column oven. Isocratic elution consisted of acetonitrile:water (FA 0.5%) (4:1) was done in 30min. The flow rate was maintained at 0.8 ml/min. The cannabinoids CBD, CBG, CBN and THC were monitored at 225 and CBDA, CBGA were monitored at 306 nm respectively using dual absorbance detector Waters 2487 (Milford, MA, USA). The injection volume of 0.1 mg/ml sample was 10 µl. Data evaluation was performed using Clarity software.

Quantification of cannabinoids was obtained from linear regression equation of calibration curve of individual reference standard by plotting concentration versus the area ratio

Analysis of terpenes was performed using GC-HS system equipped with auto injector. Separation was accomplished on Rx624Sil, 30m, 0.25 µm ID column.